


**Asilomed**<sup>®</sup>  
**PRIVATE LIMITED** 

**ISO 13485:2016 Certified Company**

**All Products FDA Approved**

**Product Guide**

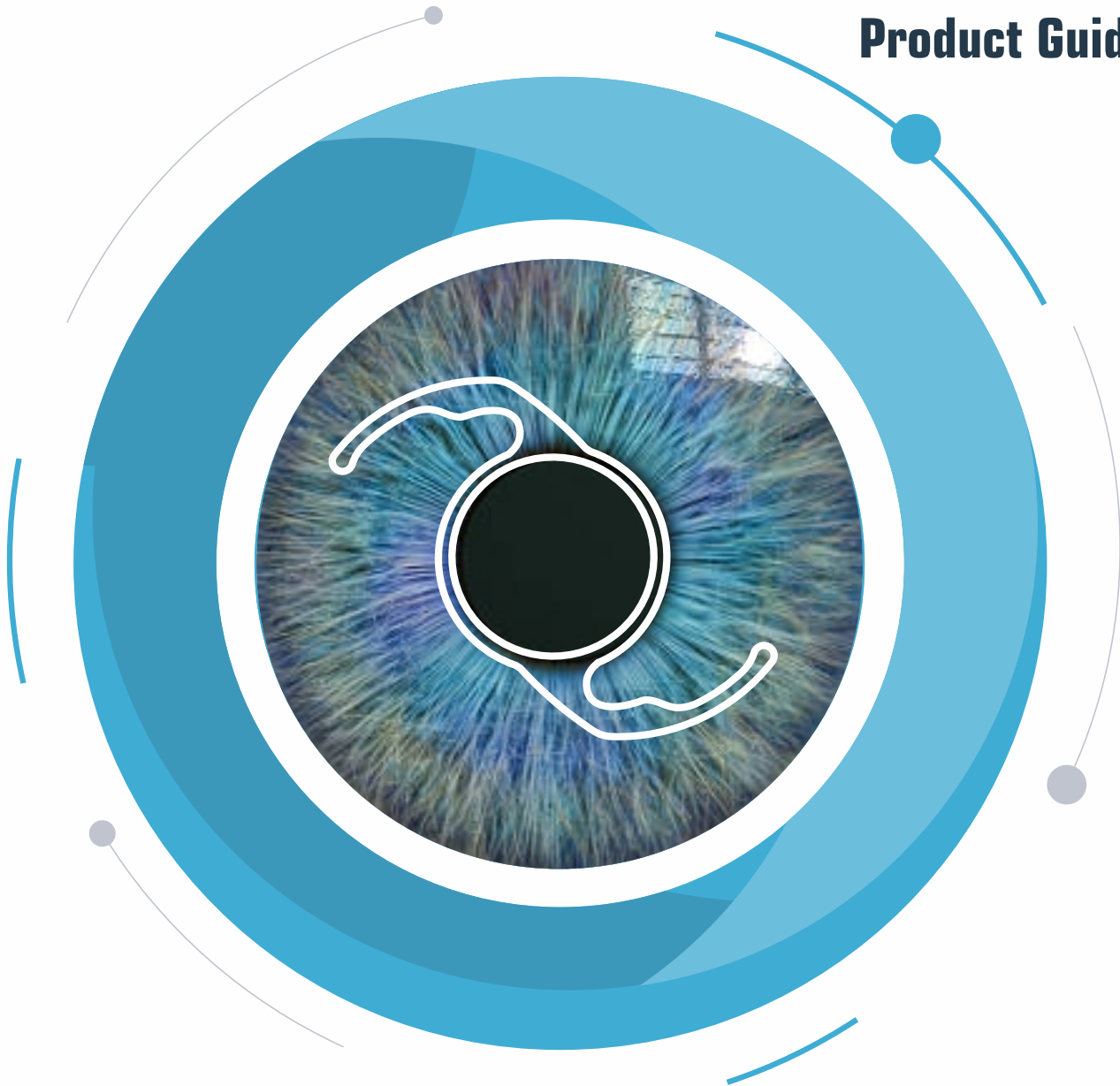

**Manufacturers of**  
All kinds of Intraocular Lens  
Ophthalmic Surgical Blades  
Ophthalmic Solutions  
Ophthalmic Drapes

**Raw Material Partner**  
**Contamac**  
EMPOWERING YOUR VISION

[www.asilomed.com](http://www.asilomed.com)

Catalogue 01-R0/03.01.2024

Welcome To



**ASILOMED PRIVATE LIMITED** is reputed Manufacturer and Exporter of Medical Devices the products include Intraocular Lenses, Ophthalmic Surgical Blades, Ophthalmic Solutions & Ophthalmic Drapes.

**ABOUT ASILOMED:**

Our Company ASILOMED was originally found in the year of 2020 by our group of directors, located in Krishnagiri District, Tamilnadu, INDIA with the ambition to create innovative and high-quality products, balanced with clear strategy, ethical Considerations and focus on long term sustainability to ensure lasting Success by our dedicated team , works round the clock to bringing most modern state of art technology to manufacture high quality premium Intraocular Lens with latest world class modern clean room facility covering 30000Sq Ft area as per INDIAN and EU MDR .

**CERTIFICATIONS:**

Asilomed is an ISO 13485:2016 certified company for Design, Development, Manufacture and Supply of All Kinds of Intraocular Lenses, Ophthalmic Surgicals Knives, Ophthalmic Solutions, and Ophthalmic Drapes. We are also a licensed manufacturer of Class A, B and Class C Medical Devices from Central Drugs Standard Control Organization (CDSCO).

**KEY SERVICES AND ASPECTS:**

We do undertake contract Manufacturing (OEM) and Other Brand Labels (OBL) for people becoming our business partners across the world to benefit competitive advantage in their respective regions and markets by the their association with us with the key aspects of

- Focus on Quality Management
- Compliance with Regulatory Requirements
- Risk Management
- Traceability and Documentation
- Continuous Improvement
- Setting Standards
- Monitoring and Surveillance



**Company profile**

Company incorporated in Krishnagiri, Tamilnadu, India in the year of 2020.

Our IOL company established as per Indian & European New Medical Device Regulations ( MDR).




# 3 Piece

## Hydrophobic Acrylic 3 Piece-IOL

### AMH 6130M

**Asilo 3 Piece** is Clearest Hydrophobic Acrylic Including the clearest vacuole-free material of all foldable IOLs tested. Its stable, **Asilo 3 Piece** architecture is designed for long-term centration and refractive stability. **Asilo 3 Piece IOL** reduces glare and its square posterior edge provides a continuous 360° capsular contact. With the **Asilo Jet™** implantation system, the **Asilo 3 Piece IOL** can be delivered through an incision as small as **2.8 mm**

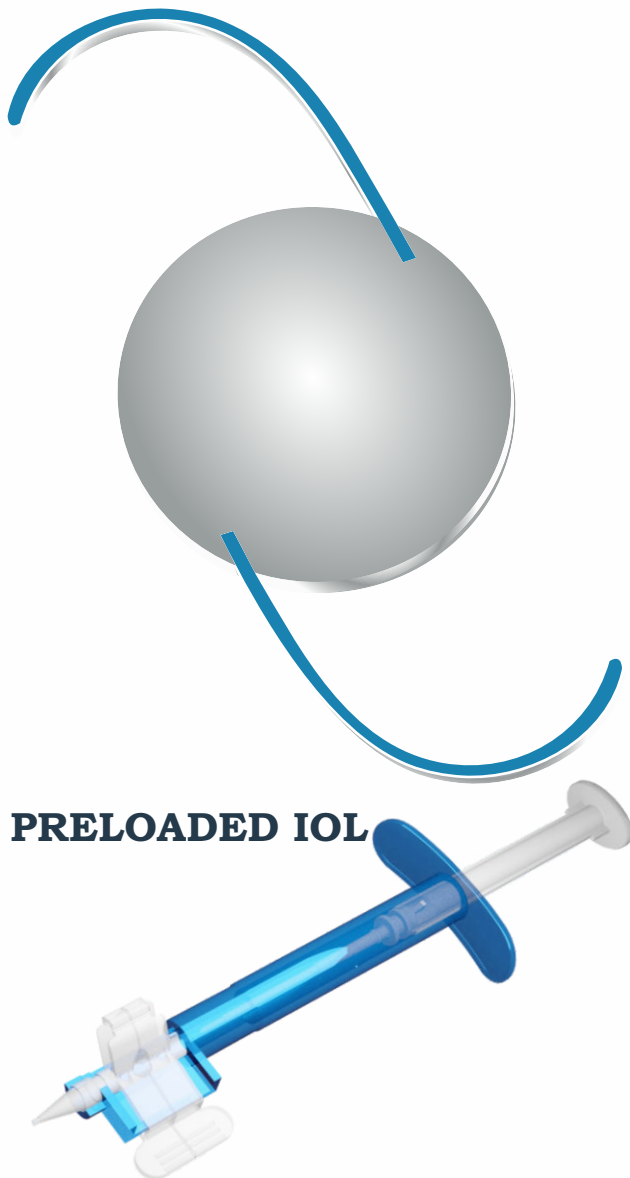

**PRELOADED IOL**

#### POWER RANGE

##### Standard Powers

10.00 D → +30.0 D (0.5 D steps)

##### Haptic Material

PVDF ( Polyvinylidene fluoride Monofilament)  
**Break free Haptics**

#### MATERIAL

Hydrophobic Acrylic blanks contain an UV Absorber with <0.5% Water Content. optional blue light filtering Chromophore

##### Refractive Index

1.493

##### Sterilization

Ethylene Oxide

#### DESIGN

##### Optic Diameter Overall Diameter

6 mm

13.00 mm, 13.50 mm

##### Optic Design

Biconvex - 360° Square Edge with Aspheric Optics

##### Nominal A-constant

118.7

##### Haptic Design

Modified "C" Loop, 0° Angulation

##### PCO Prevention

360° Square Edge

##### Shelf Life

5 years after sterilization

#### INJECTION

##### Estimated Incision Size

2.6 - 2.8 mm

##### Injector & Cartridge

**Asilo Jet™** -Front Load with 2.8 mm Cartridge

#### CONSTANTS

*for optical and immersion ultrasound measurements*

##### SRK/T A-Constant

A = 118.7\*

##### Haigis

a0 = 1.260\*  
a1 = 0.400\*  
a2 = 0.100\*

##### Hoffer Q

pACD = 5.39\*

##### Holladay I

SF = 1.85\*

##### SRK/ II A-Constant

A = 118.9\*

**\*Optimized constants**



# 3 Piece

## PMMA 3 Piece IOL

### MP 6125H

**Asilo 3 Piece IOL** optic with an accuracy in the nanometer range lays the foundation for premium quality and produces brilliant, clear, and sharp images. A high level of stability and centering ensures the brilliant imaging properties of the optic.

The modified C-loop haptics made of highly flexible, blue-colored **PVDF ( Polyvinylidene fluoride Monofilament)** ensure optimal IOL centering thanks to high axial stability and radial flexibility. In addition, the fine haptics allow the anterior and posterior capsule membranes to merge for additional stability. The fine C-loop haptic is also suitable for the sensitive sulcus tissue and minimizes irritation due to the small contact area. The optic material made of clinically proven acrylat is characterized by excellent imaging quality. It is completely glistening-free and demonstrates excellent uveal biocompatibility.

#### POWER RANGE

##### Standard Powers

10.00 D → +30.0 D (0.5 D steps)

##### Haptic Material

PVDF ( Polyvinylidene fluoride Monofilament)

#### MATERIAL

PMMA 3 Piece Intraocular Lens with break free Haptics

##### Refractive Index

1.495

##### Sterilization

Ethylene Oxide

#### DESIGN

##### Optic Diameter Overall Diameter

6 mm

13.00 mm, 13.50 mm

Biconvex - 360° Square Edge with Aspheric Optics

##### Nominal A-constant

118.4

##### Haptic Design


Modified "C" Loop, 5° Angulation

##### PCO Prevention

360° Square Edge

##### Shelf Life

5 years after sterilization

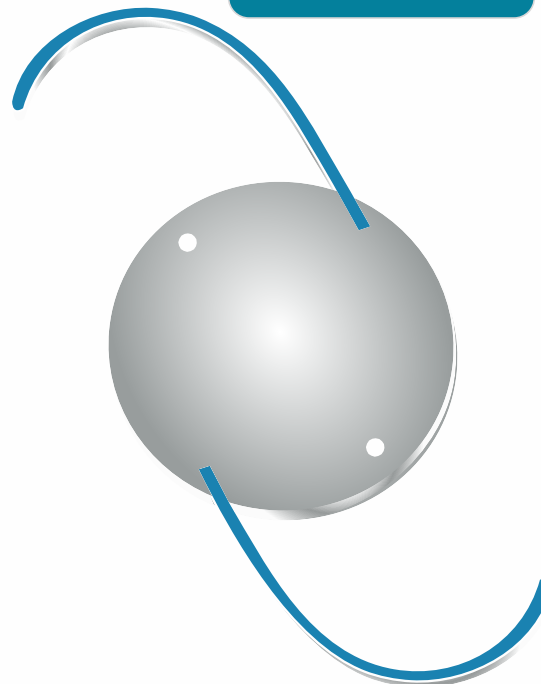



# ED-FOCUS™


## Hydrophobic Acrylic EDOF-IOL

### AMH-EDF-6125

The extended depth of focus IOL ( **EDOF IOL**) IOL is a single, contiguous, elongated focal point that enhances depth of focus

The extended depth of focus IOL ( EDOF IOL) Extended Range of Vision, is a new technology that has recently emerged in the treatment of Presbyopia-correcting IOLs.

It is based on the principle of creating a single elongated focal point to enhance the “range of vision” or “depth of focus”.

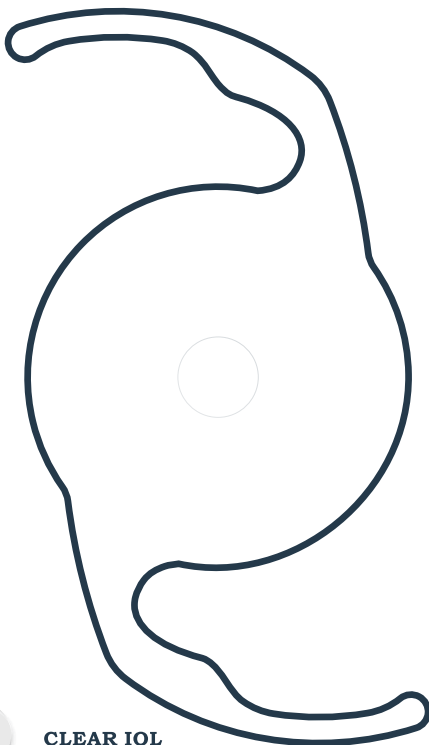

 CLEAR IOL

 YELLOW IOL

**PRELOADED IOL**



#### POWER RANGE

##### Standard Powers

5.00 D → +30.0 D (0.5 D steps)

##### Intermediate Addition

2.5 D

#### MATERIAL

Hydrophobic Acrylic blanks contain an UV Absorber with <0.5% Water Content. optional blue light filtering Chromophore

##### Refractive Index

1.493

##### Sterilization

Ethylene Oxide

#### DESIGN

##### Optic Diameter Overall Diameter

6 mm

13.00 mm

##### Optic Design

Biconvex - 360° Square Edge with Refractive Optics

##### Nominal A-constant

118.7

##### Haptic Design

Single Haptics, 0° Angulation

##### PCO Prevention

360° Square Edge

##### Shelf Life

5 years after sterilization

#### INJECTION

##### Estimated Incision Size

2.4 - 2.6 mm

##### Injector & Cartridge

**Asilo Jet™** -Front Load with 2.4 mm Cartridge

#### CONSTANTS

*for optical and immersion ultrasound measurements*

##### SRK/T A-Constant

A = 118.7\*

##### Haigis

a0 = 1.260\*  
a1 = 0.400\*  
a2 = 0.100\*

##### Hoffer Q

pACD = 5.39\*

##### Holladay I

SF = 1.85\*

##### SRK/ II A-Constant

A = 118.9\*

\*Optimized constants



**TRIO<sup>®</sup>**


**Hydrophobic Acrylic Trifocal -IOL**

**AMH-TRI-6125**

Single Piece Sterile with 360° Square Edge **TRIFOCAL** Hydrophobic Acrylic Foldable Aspheric Posterior Chamber Intraocular Lens(IOL)with UV blocker, Optionally with blue light filter.

This optic is **TRIFOCAL**

The IOL is non-Preloaded , which has to be manually loaded into a compatible Injector.

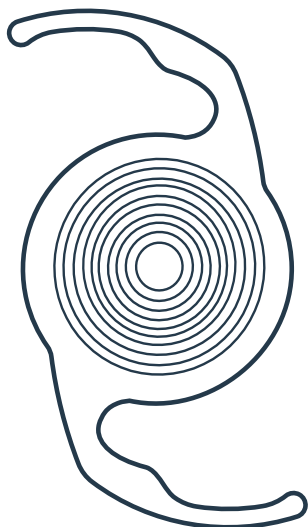

 CLEAR IOL

 YELLOW IOL

**PRELOADED IOL**



**POWER RANGE**

**Standard Powers** 0.00 D → +35.0 D (0.5 D steps)

**Intermediate Addition** +3.50 D Near  
+1.75 D Intermediate

**MATERIAL**

Hydrophobic Acrylic blanks contain an UV Absorber with <0.5% Water Content, optional blue light filtering Chromophore

**Refractive Index** 1.493

**Sterilization** Ethylene Oxide

**DESIGN**

**Optic Diameter Overall Diameter** 6 mm 12.50 mm  
13.00 mm

**Optic Design** Aspheric - Aberration Neutral  
Biconvex (+0.00 D to 35.00 D)

**Nominal A-constant** 118.7

**Haptic Design** Single Haptics with 0° Angulation

**PCO Prevention** 360° Square Edge

**Shelf Life** 5 years after sterilization

**INJECTION**

**Estimated Incision Size** 2.4 - 2.6 mm

**Injector & Cartridge** **Asilo Jet™** -Front Load with 2.4 mm Cartridge

**CONSTANTS**  
*for optical and immersion ultrasound measurements*

SRK/T A-Constant	Haigis	Hoffer Q
A = 118.7*	a0 = 1.260* a1 = 0.400* a2 = 0.100*	pACD = 5.39*

Holladay I	SRK/ II A-Constant
SF = 1.85*	A = 118.9*

**\*Optimized constants**




## Hydrophobic Acrylic Foldable -IOL

**AMH 6125SQ/  
AMH 6130SQ**

Single Piece Sterile with 360° Square Edge  
**HYDROPHOBIC ASPHERIC** Acrylic Foldable  
Aspheric Posterior Chamber Intraocular  
Lens(IOL)with UV blocker, Optionally with  
blue light filter.

This optic is **MONOFOCAL**

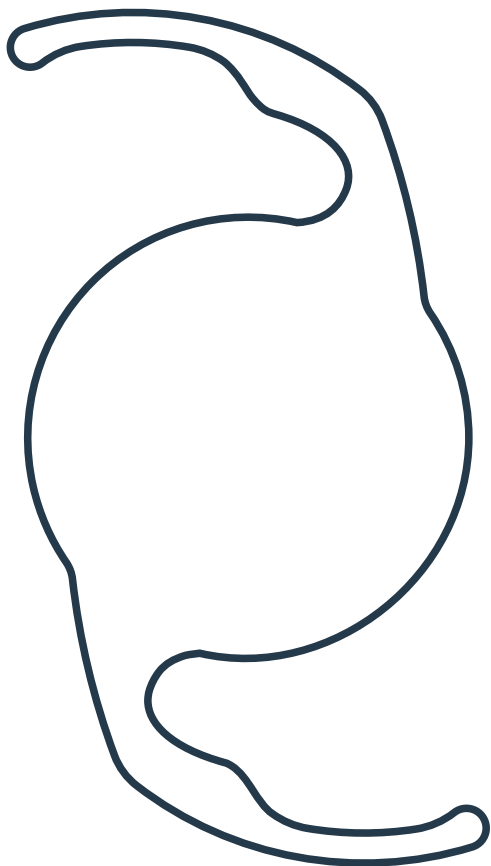

**PRELOADED IOL**



### POWER RANGE

#### Standard Powers

00.0 D → +30.0 D (0.5 D steps)

#### Cylinders

-10.0 D → -1.0 D (0.5 D steps)  
+31.0 D → +45.0 D (0.5 D steps)

### MATERIAL

Hydrophobic Acrylic blanks contain an UV Absorber with <0.5% Water Content,  
optional blue light filtering Chromophore

#### Refractive Index

1.493

#### Sterilization

Ethylene Oxide

### DESIGN

#### Optic Diameter Overall Diameter

6 mm

12.50 mm & 13.00 mm

#### Optic Design

Aspheric - Aberration Neutral  
Biconvex (+0.00 D to 35.00 D)

#### Nominal A-constant

118.7

#### Haptic Design

Single Haptics with 0° Angulation

#### PCO Prevention

360° Square Edge

#### Shelf Life

5 years after sterilization

### INJECTION

#### Estimated Incision Size

2.4 - 2.6 mm

#### Injector & Cartridge

Front Load with 2.4 mm Cartridge

### CONSTANTS

*for optical and immersion ultrasound measurements*

#### SRK/T A-Constant

A = 118.7\*

#### Haigis

a0 = 1.260\*  
a1 = 0.400\*  
a2 = 0.100\*

#### Hoffer Q

pACD = 5.39\*

#### Holladay I

SF = 1.85\*

#### SRK/ II A-Constant

A = 118.9\*

\*Optimized constants



**VUE+**<sup>®</sup>


**YELLOW Hydrophobic Acrylic Foldable -IOL**

**AMH 6125SQY /  
AMH 6130SQY**

Single Piece Sterile with 360° Square Edge  
**HYDROPHOBIC ASPHERIC** Acrylic Foldable  
Aspheric Posterior Chamber Intraocular  
Lens(IOL) with UV blocker, Optionally with  
blue light filter.

This optic is **MONOFOCAL**

The IOL is non-Preloaded , which has  
to be manually loaded into a compatible  
Injector.

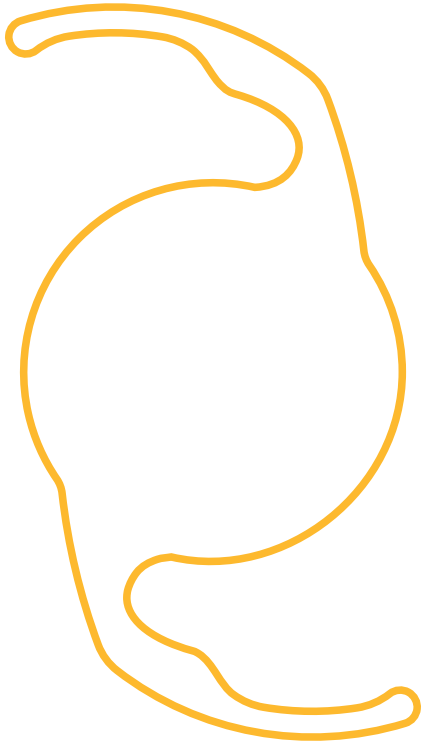

**PRELOADED IOL**



**POWER RANGE**

<b>Standard Powers</b>	00.0 D → +30.0 D (0.5 D steps)
<b>Extreme Powers</b>	-10.0 D → -1.0 D (0.5 D steps) +31.0 D → +45.0 D (0.5 D steps)

**MATERIAL**

Hydrophobic Acrylic blanks contain an UV Absorber with <0.5% Water Content.  
optional blue light filtering Chromophore

<b>Refractive Index</b>	1.493
<b>Sterilization</b>	Ethylene Oxide

**DESIGN**

<b>Optic Diameter Overall Diameter</b>	6 mm	13.00 mm
<b>Optic Design</b>	Aspheric - Aberration Neutral Biconvex (+0.00 D to 35.00 D)	
<b>Nominal A-constant</b>	118.7	
<b>Haptic Design</b>	Single Haptics with 0° Angulation	
<b>PCO Prevention</b>	360° Square Edge	
<b>Shelf Life</b>	5 years after sterilization	

**INJECTION**

<b>Estimated Incision Size</b>	2.4 - 2.6 mm
<b>Injector &amp; Cartridge</b>	<b>Asilo Jet™</b> -Front Load with 2.4 mm Cartridge

**CONSTANTS**  
*for optical and immersion ultrasound measurements*

SRK/T A-Constant	Haigis	Hoffer Q
A = 118.7*	a0 = 1.260* a1 = 0.400* a2 = 0.100*	pACD = 5.39*
Holladay I	SRK/ II A-Constant	
SF = 1.85*	A = 118.9*	

**\*Optimized constants**





# FOLD<sup>®</sup>

## Hydrophilic Acrylic Foldable-IOL

Single Piece Sterile with 360° Square Edge

Hydrophilic Acrylic Foldable



**SPHERICAL** Posterior Chamber

Intraocular Lens(IOL) with UV blocker,

Optionally with blue light filter.

This optic is **MONOFOCAL**

The IOL is non-Preloaded , which has to be manually loaded into a compatible Injector.

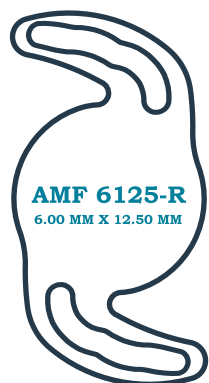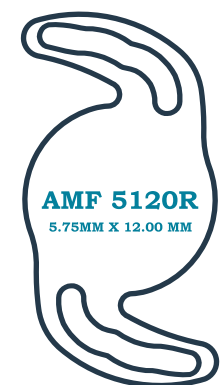

### POWER RANGE

#### Standard Powers

00.0 D → +30.0 D (0.5 D steps)

#### Extreme Powers

-10.0 D → -1.0 D (0.5 D steps)  
+31.0 D → +45.0 D (0.5 D steps)

### MATERIAL

Copolymer of hydroxyethyl methacrylate and methyl methacrylate with 26% water content, UV blocker and optional blue light filter

#### Refractive Index

1.46

#### Sterilization

Steam

### DESIGN

#### Optic Design

Negative Aspheric - Aberration Neutral  
Convex-concave (-10.0 D → -1.0 D)  
Biconvex (0.0 D → +45.0 D)

#### Nominal A-constant

118.0

#### Haptic Design

Double Haptics, 0° Angulation

#### PCO Prevention

360° Square Edge

#### Shelf Life

3 years after sterilization

### INJECTION

#### Estimated Incision Size

2.2 - 2.4 mm

#### Injector & Cartridge

**Asilo Jet™** -Front Load with 2.2 mm Cartridge

### CONSTANTS

*for optical and immersion ultrasound measurements*

#### SRK/T A-Constant

A = 118.0\*

#### Haigis

a0 = 0.568\*  
a1 = 0.400\*  
a2 = 0.100\*

#### Hoffer Q

pACD = 4.80\*

#### Holladay I

SF = 1.02\*

#### SRK/ II A-Constant

A = 117.88

\*Optimized constants



**FLEX<sup>®</sup>**

## Hydrophilic Aspheric Foldable-IOL


Single Piece Sterile with 360° Square Edge  
Hydrophilic Acrylic Foldable

**CLEAR ASPHERIC** Posterior Chamber

Intraocular Lens(IOL)with UV  
blocker, Optionally with blue light filter.

This optic is **MONOFOCAL**

The IOL is non-Preloaded , which has  
to be manually loaded into a compatible  
Injector.

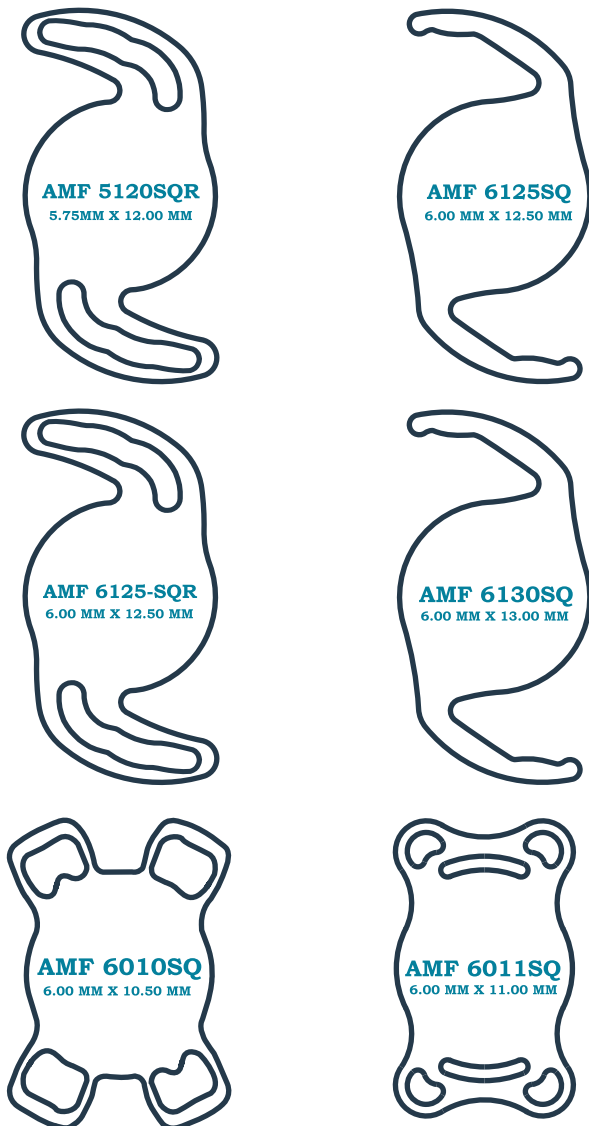

POWER RANGE		
<b>Standard Powers</b>	00.0 D → +30.0 D (0.5 D steps)	
<b>Extreme Powers</b>	-10.0 D → -1.0 D (0.5 D steps) +31.0 D → +45.0 D (0.5 D steps)	
MATERIAL		
Copolymer of hydroxyethyl methacrylate and methyl methacrylate with 26% water content, UV blocker and optional blue light filter		
<b>Refractive Index</b>	1.46	
<b>Sterilization</b>	Steam	
DESIGN		
<b>Optic Design</b>	Negative Aspheric - Aberration Neutral Convex-concave (-10.0 D → -1.0 D) Biconvex (0.0 D → +45.0 D)	
<b>Nominal A-constant</b>	118.0	
<b>Haptic Design</b>	Single Haptics, 0° Angulation	
<b>PCO Prevention</b>	360° Square Edge	
<b>Shelf Life</b>	3 years after sterilization	
INJECTION		
<b>Estimated Incision Size</b>	2.2 - 2.4 mm	
<b>Injector &amp; Cartridge</b>	Asilo Jet™ -Front Load with 2.2 mm Cartridge	
CONSTANTS <i>for optical and immersion ultrasound measurements</i>		
<b>SRK/T A-Constant</b>	<b>Haigis</b>	<b>Hoffer Q</b>
A = 118.0*	a0 = 0.568* a1 = 0.400* a2 = 0.100*	pACD = 4.80*
<b>Holladay I</b>	<b>SRK/ II A-Constant</b>	
SF = 1.02*	A = 117.88	

\*Optimized constants



# TAN™

## Hydrophilic Yellow Aspheric Foldable-IOL

Single Piece Sterile with 360° Square Edge  
 Hydrophilic Acrylic Foldable  
**YELLOW ASPHERIC** Posterior Chamber  
 Intraocular Lens(IOL)with UV  
 blocker, Optionally with blue light filter.

This optic is **MONOFOCAL**

The IOL is non-Preloaded , which has  
 to be manually loaded into a compatible  
 Injector.



### POWER RANGE

<b>Standard Powers</b>	00.0 D → +30.0 D (0.5 D steps)
------------------------	--------------------------------

<b>Extreme Powers</b>	-10.0 D → -1.0 D (0.5 D steps) +31.0 D → +45.0 D (0.5 D steps)
-----------------------	---

### MATERIAL

Copolymer of hydroxyethyl methacrylate and methyl methacrylate with 26% water content, UV blocker and optional blue light filter

<b>Refractive Index</b>	1.46
-------------------------	------

<b>Sterilization</b>	Steam
----------------------	-------

### DESIGN

<b>Optic Design</b>	Negative Aspheric - Aberration Neutral Convex-concave (-10.0 D → -1.0 D) Biconvex (0.0 D → +45.0 D)
---------------------	---

<b>Nominal A-constant</b>	118.0
---------------------------	-------

<b>Haptic Design</b>	Single Haptics, 0° Angulation
----------------------	-------------------------------

<b>PCO Prevention</b>	360° Square Edge
-----------------------	------------------

<b>Shelf Life</b>	3 years after sterilization
-------------------	-----------------------------

### INJECTION

<b>Estimated Incision Size</b>	2.2 - 2.4 mm
--------------------------------	--------------

<b>Injector &amp; Cartridge</b>	<b>Asilo Jet™</b> -Front Load with 2.2 mm Cartridge
---------------------------------	---

### CONSTANTS for optical and immersion ultrasound measurements

SRK/T A-Constant	Haigis	Hoffer Q
A = 118.0*	a0 = 0.568* a1 = 0.400* a2 = 0.100*	pACD = 4.80*

Holladay I	SRK/ II A-Constant
SF = 1.02*	A = 117.88

\*Optimized constants

# Asilo LENS™

## PMMA Intraocular Lens



**AMP 5125**



**AMP 5212**




**AMP 5525**



**AMP 6125H**

<b>Standard Powers</b>	00.0 D → +35.0 D (0.5 D steps)	00.0 D → +35.0 D (0.5 D steps)	00.0 D → +35.0 D (0.5 D steps)	00.0 D → +35.0 D (0.5 D steps)
<b>Extreme Powers</b>	10.0 D → 1.0 D (1.0 D steps) +35.0 D → +45.0 D (1.0 D steps)	10.0 D → 1.0 D (1.0 D steps) +35.0 D → +45.0 D (1.0 D steps)	10.0 D → 1.0 D (1.0 D steps) +35.0 D → +45.0 D (1.0 D steps)	10.0 D → 1.0 D (1.0 D steps) +35.0 D → +45.0 D (1.0 D steps)
<b>Optic Diameter</b>	5.00 mm	5.25 mm	5.50 mm	6.00 mm
<b>Overall Diameter</b>	12.50 mm	12.00 mm	12.50 mm	12.50 mm
<b>Optic Design</b>	Equiconvex	Equiconvex	Equiconvex	Equiconvex
<b>Haptic Design</b>	“Cap.C” Design with StepVault	“Mod.C” Design with StepVault	“Mod.C” Design with StepVault	“Mod.C” Design with StepVault
<b>Ac Depth</b>	5.00 mm	5.00 mm	5.00 mm	4.80 mm
<b>Dialing Holes</b>	Nil	Nil	Nil	2
<b>Suggested A Constant</b>	118.4*	118.4*	118.4*	118.2*
<b>Shelf Life</b>	5 years after sterilization	5 years after sterilization	5 years after sterilization	5 years after sterilization

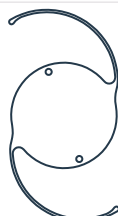

**AMP 6130H**



**AMP 6513 H**



**AMS 6125**



**ASF 65135**


<b>Standard Powers</b>	00.0 D → +35.0 D (0.5 D steps)	00.0 D → +35.0 D (0.5 D steps)	00.0 D → +35.0 D (0.5 D steps)	00.0 D → +35.0 D (0.5 D steps)
<b>Extreme Powers</b>	10.0 D → 1.0 D (1.0 D steps) +35.0 D → +45.0 D (1.0 D steps)	10.0 D → 1.0 D (1.0 D steps) +35.0 D → +45.0 D (1.0 D steps)	10.0 D → 1.0 D (1.0 D steps) +35.0 D → +45.0 D (1.0 D steps)	10.0 D → 1.0 D (1.0 D steps) +35.0 D → +45.0 D (1.0 D steps)
<b>Optic Diameter</b>	6.00 mm	6.50 mm	6.00 mm	6.50 mm
<b>Overall Diameter</b>	13.00 mm	13.00 mm	12.50 mm	13.50 mm
<b>Optic Design</b>	Equiconvex	Equiconvex	Spherical Biconvex	Spherical Biconvex
<b>Haptic Design</b>	“Mod.C” Design with Step Vault	“Mod.C” Design with Step Vault	“Ant.S” Design with Step Vault	“Mod.C” Design with Step Vault
<b>Ac Depth</b>	4.80 mm	4.80 mm	2.95 mm	4.80 mm
<b>Dialing Holes</b>	2	2	Nil	2
<b>Suggested A Constant</b>	118.2*	118.2*	115.3*	118.2*
<b>Shelf Life</b>	5 years after sterilization	5 years after sterilization	5 years after sterilization	5 years after sterilization

# Asilo Edge<sup>®</sup>


## PMMA Aspheric Optics with 360° Square Edge Intraocular Lens



**AMP 5125SQ**

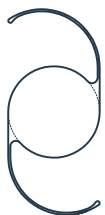

**AMP 5525SQ**

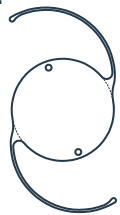


**AMP 6125HSQ**




**AMP 6130HSQ**

<b>Standard Powers</b>	00.0 D → +35.0 D (0.5 D steps)	00.0 D → +35.0 D (0.5 D steps)	00.0 D → +35.0 D (0.5 D steps)	00.0 D → +35.0 D (0.5 D steps)
<b>Extreme Powers</b>	10.0 D → 1.0 D (1.0 D steps) +35.0 D → +45.0 D (1.0 D steps)	10.0 D → 1.0 D (1.0 D steps) +35.0 D → +45.0 D (1.0 D steps)	10.0 D → 1.0 D (1.0 D steps) +35.0 D → +45.0 D (1.0 D steps)	10.0 D → 1.0 D (1.0 D steps) +35.0 D → +45.0 D (1.0 D steps)
<b>Optic Diameter</b>	5.00 mm	5.50 mm	6.00 mm	6.00 mm
<b>Overall Diameter</b>	12.50 mm	12.50 mm	12.50 mm	13.00 mm
<b>Optic Design</b>	Aberration Neutral, Negative Aspherical Biconvex	Aberration Neutral, Negative Aspherical Biconvex	Aberration Neutral, Negative Aspherical Biconvex	Aberration Neutral, Negative Aspherical Biconvex
<b>Haptic Design</b>	“Mod.C” Design with StepVault	“Mod.C” Design with StepVault	“Mod.C” Design with StepVault	“Mod.C” Design with StepVault
<b>Ac Depth</b>	5.00 mm	5.00 mm	4.80 mm	4.80 mm
<b>Dialing Holes</b>	Nil	Nil	2	2
<b>Suggested A Constant</b>	118.4*	118.4*	118.4*	118.4*
<b>Shelf Life</b>	5 years after sterilization	5 years after sterilization	5 years after sterilization	5 years after sterilization

# Asilo IRIS<sup>®</sup> PMMA IRIS CLAW Intraocular Lens

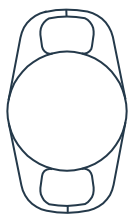

**AIC 4272**

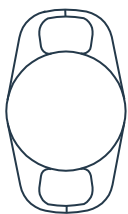

**AIC 5085**



**AIC 5580**



**AIC 5590**

<b>Standard Powers</b>	00.0 D → +35.0 D (0.5 D steps)	00.0 D → +35.0 D (0.5 D steps)	00.0 D → +35.0 D (0.5 D steps)	00.0 D → +35.0 D (0.5 D steps)
<b>Optic Diameter</b>	4.00 mm	5.00 mm	5.50 mm	5.50 mm
<b>Overall Diameter</b>	7.25 mm	8.50 mm	8.50 mm	9.00 mm
<b>Optic Design</b>	Biconvex Convex Concave (according to Diopter)	Biconvex Convex Concave (according to Diopter)	Biconvex Convex Concave (according to Diopter)	Biconvex Convex Concave (according to Diopter)
<b>Haptic Design</b>	“Mod.C” Design with StepVault	“Mod.C” Design with StepVault	“Mod.C” Design with StepVault	“Mod.C” Design with StepVault
<b>Ac Depth</b>	2.95 mm	2.95 mm	2.95 mm	2.95 mm
<b>Suggested A Constant</b>	Anterior Side- 115.0* Posterior Side- 117.4*	Anterior Side- 115.0* Posterior Side- 117.4*	Anterior Side- 115.0* Posterior Side- 117.4*	Anterior Side- 115.0* Posterior Side- 117.4*
<b>Shelf Life</b>	5 years after sterilization	5 years after sterilization	5 years after sterilization	5 years after sterilization




**Gold™**


**YELLOW PMMA Aspheric Optics with 360° Square Edge Intraocular Lens**




**AMP 5125SQY**

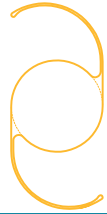

**AMP 5525SQY**

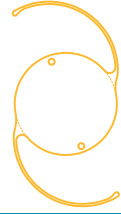

**AMP 6125HSQY**

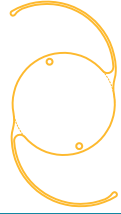


**AMP 6130HSQY**

<b>Standard Powers</b>	00.0 D → +35.0 D (0.5 D steps)	00.0 D → +35.0 D (0.5 D steps)	00.0 D → +35.0 D (0.5 D steps)	00.0 D → +35.0 D (0.5 D steps)
<b>Extreme Powers</b>	10.0 D → 1.0 D (1.0 D steps) +35.0 D → +45.0 D (1.0 D steps)	10.0 D → 1.0 D (1.0 D steps) +35.0 D → +45.0 D (1.0 D steps)	10.0 D → 1.0 D (1.0 D steps) +35.0 D → +45.0 D (1.0 D steps)	10.0 D → 1.0 D (1.0 D steps) +35.0 D → +45.0 D (1.0 D steps)
<b>Optic Diameter</b>	5.00 mm	5.50 mm	6.00 mm	6.00 mm
<b>Overall Diameter</b>	12.50 mm	12.50 mm	12.50 mm	13.00 mm
<b>Optic Design</b>	Aberration Neutral, Negative Aspherical Biconvex	Aberration Neutral, Negative Aspherical Biconvex	Aberration Neutral, Negative Aspherical Biconvex	Aberration Neutral, Negative Aspherical Biconvex
<b>Haptic Design</b>	“Mod.C” Design with StepVault	“Mod.C” Design with StepVault	“Mod.C” Design with StepVault	“Mod.C” Design with StepVault
<b>Ac Depth</b>	5.00 mm	5.00 mm	4.80 mm	4.80 mm
<b>Dialing Holes</b>	Nil	Nil	2	2
<b>Suggested A Constant</b>	118.4*	118.4*	118.4*	118.4*
<b>Shelf Life</b>	5 years after sterilization	5 years after sterilization	5 years after sterilization	5 years after sterilization




**RING™**


**PMMA Capsular Tension Ring**

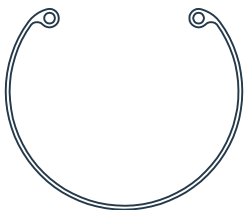


**AMR 1109**

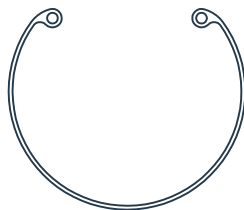

**AMR 1210**

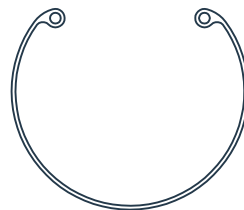

**AMR 1311**

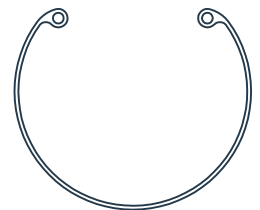

**AMR 1412**

<b>Overall Diameter</b>	11.00 mm	12.00 mm	13.00 mm	14.00 mm
<b>Insertion Diameter</b>	9.00 mm	10.00 mm	11.00 mm	12.00 mm
<b>Material</b>	CQ-PMMA Clear/Blue	CQ-PMMA Clear/Blue	CQ-PMMA Clear/Blue	CQ-PMMA Clear/Blue



**AMR 1109L / AMR 1210 L  
AMR 1311L / AMR 1412 L**



**AMR 1109R / AMR 1210 R  
AMR 1311R / AMR 1412 R**



**AMR 1109LR / AMR 1210 LR**





**AMR 1311LR / AMR 1412 LR**


<b>Overall Diameter</b>	11.00 mm/12.00 mm / 13.00 mm/14.00 mm	11.00 mm/12.00 mm / 13.00 mm/14.00 mm	11.00 mm/12.00 mm	13.00 mm/14.00 mm
<b>Insertion Diameter</b>	09.00 mm/10.00 mm / 11.00 mm/12.00 mm	09.00 mm/10.00 mm / 11.00 mm/12.00 mm	09.00 mm/10.00 mm	11.00 mm/12.00 mm
<b>Fixation Hook</b>	Left side with Step Vault	Right side with Step Vault	Both side with Step Vault	Both side with Step Vault
<b>Material</b>	CQ-PMMA Clear/Blue	CQ-PMMA Clear/Blue	CQ-PMMA Clear/Blue	CQ-PMMA Clear/Blue


# Asilo <sup>TM</sup> Glide


## Ophthalmic Micro Surgical Blades

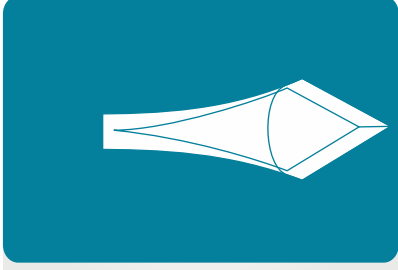
KERATOME (Sharp Tip)	MODEL	SIZE	DESCRIPTION
	AG2415	1.50mm to 1.75mm	24 Guage
	AG2018	1.80mm to 1.95mm	20 Guage
	AG2020	2.00mm to 2.25mm	20 Guage
	AG2030	2.30mm to 2.55mm	20 Guage
	AG2060	2.60mm to 2.75mm	20 Guage
	AG2028	2.80mm to 2.95mm	20 Guage
	AG2430	3.00mm to 3.20mm	24 Guage

CRESCENT	MODEL	SIZE	DESCRIPTION
	AG2020 BU	2.00mm to 2.25mm	20 Guage
	AG2020 BD	2.00mm to 2.25mm	20 Guage

KERATOME (Blunt Tip)	MODEL	SIZE	DESCRIPTION
	AG1835	3.50mm to 3.75mm	18 Guage
	AG1651	5.10mm to 5.25mm	16 Guage
	AG1655	5.50mm to 5.75mm	16 Guage

LANCE TIP	MODEL	SIZE	DESCRIPTION
	AG2415	15°±1°	24 Guage
	AG2430	30°±1°	24 Guage

MVR BLADES	MODEL	SIZE	DESCRIPTION
	AG1996	1.35mm±0.05 mm	19 Guage
	AG2096	2.00mm±0.05 mm	20 Guage
	AG2496	0.85mm±0.05 mm	24 Guage

SEGMENT RING	MODEL	AMR - 09
	Insertion Diameter	09.00 mm
	Eyelets	3
	Fixation Hook	0.40 mm
	Material	CQ-PMMA Clear/Blue

# Asilo Aqua™

Hydroxy Propyl Methyl Cellulose Ophthalmic Solution USP 2% w/v



# Asilo Hyal™

Sodium Hyaluronate Ophthalmic Solution with Prefilled Syringe





# Asilo Blue™

Trypan Blue Ophthalmic Solution



# Asilo Chol™

Carbachol Intraocular Solution USP

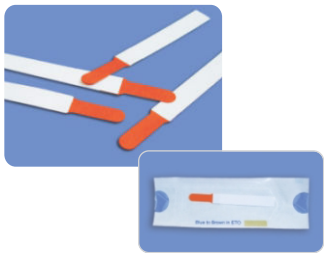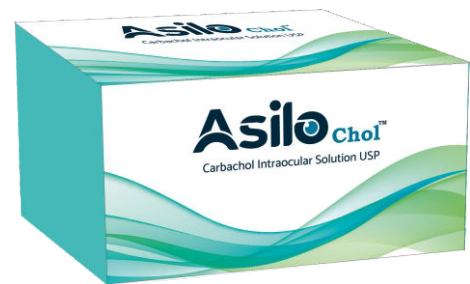

## Fluorescein Strips


- ❖ Useful for evaluating hard contact lens fitting
- ❖ Evaluating tear film breakup time
- ❖ To visualise defects in the corneal epithelium
- ❖ Each strip is impregnated with approx. 1.0 mg of Fluorescein sodium

Asilo FLURO™  
Fluorescein Sodium Ophthalmic Strips USP



## Schirmer Strips

- ❖ Used for measurement of tear production
- ❖ Each strip has a printed mm scale so that they can be easily read
- ❖ Innovative peel open packing with transparent film on one side

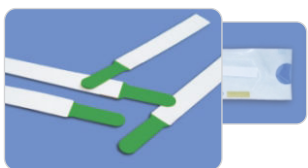

## Lissamine Green Strips

- ❖ Used to stain preferentially damaged or devitalized cells
- ❖ To indicate any dry patches, as well as any mucus-deficient or damaged corneal epithelial cells
- ❖ Each strip is impregnated with approximately 1.50 mg of Lissamine Green

## Rose Bengal Strips

- ❖ It stains devitalized or degenerated epithelial cells as well as mucous filaments
- ❖ Each strip is impregnated with approximately 1.50 mg of Rose Bengal
- ❖ Innovative peel open packing with transparent film on one side







### SMS-EYE DRAPE WITH FLUID COLLECTION BAG

S.No	CODE	DESCRIPTION
1	AD 5040	50 X 40 CM
2	AD 6060	60 X 60 CM
3	AD 7070	70 X 70 CM
4	AD 7080	70 X 80 CM
5	AD 8090	80 X 90 CM
6	AD 100100	100 X 100 CM
7	AD 100120	100 X 120 CM
8	AD 120150	120 X 150 CM
9	AD 150150	150 X 150 CM

### CONSUMABLES - DISPOSABLE

S.No	CODE	DESCRIPTION
1	ASG-001	SURGEON GOWN
1	APG-002	PATIENT GOWN
1	AEP-003	EYE PAD
1	APK-004	PERSONAL PROTECTIVE KIT
1	ASC-005	SURGEON CAP
1	ASC-006	SHOE COVER
1	ABP-007	BED SHEET WITH PILLOW COVER
1	TBS-008	BED SHEET - SMALL
1	ABS-009	BED SHEET - BIG
1	AEP-010	EYE PAD WITH STICKER

### OPHTHALMIC SURGERY KIT

S.No	CONTENTS	QTY
1	Eye Drape - 100 x 100 CM	1No
2	Trolley Drapes	1No
3	Eye Pad With Sticker	1No
4	Speculum	1No
5	Irrigating Cystotome	1No
6	Hydrodissection Cannula	1No
7	Disposable Simcoe Cannula	1No
8	J & J Buds	1No
9	Swab for Cotton Balls	1No
10	Syringes	1No
11	Air Cannula	1No
12	26 G Needles	1No
13	SICS Blades	1No
14	Surgeon Gown	1No
15	Patient Gown	1No
16	Shoe Cover	1No
17	PATIENT BED COVER	1No
18	NASAL TUBES	1No
19	Plastic Bowels	1No

### LSCS KIT

S.No	CONTENTS	QTY
1	Surgeon Gown	3No
2	Hand Towel	3No
3	Centre Hole Drape	1No
4	Medium Sheet	4No
5	Large Sheet	2No
6	Wrapper	1No
7	Apron Large	1No
8	Carry Bag	1No
9	Baby Towel	1No
10	Trolley Sheet	1No

### LSCS KIT

S.No	CONTENTS	QTY
1	Patient Gown	1No
2	Patient Cap	1No
3	Leggings	1Pair

**Get in touch now**

**Mfg.Lic No: MFG/MD/2023/000633**

**Manufacturing Unit:**

Survey No.209/1B, Jinjupalli Village,  
Billanakuppam Post, Krishnagiri-635115,  
Tamilnadu , INDIA. Tel: +91 4343 -291770  
Mob: +91-9788852516 / + 91 - 8220049863  
Email- admin@asilomed.com , info@asilomed.com



EC	REP
----	-----

Europecert (UG) (haftungsbeschränkt)  
Alsstr. 97, 41063 Mönchengladbach,  
Germany

 **Technology By:**

**OCUSERV INSTRUMENTS INC.,**  
No.1459, East 27th Street,  
NY BROOKLYN- 11210,  
U.S.A.

Email- [techd@asilomed.com](mailto:techd@asilomed.com)



Website: [www.asilomed.com](http://www.asilomed.com)